How do I inject with the Pre-filled Pen?

Yuflyma[®] is offered in an easy-to-use pre-filled pen device. The following steps will guide you to a safe and easy administration of Yuflyma®.

Important information

- Use the pre-filled pen only if your healthcare provider has trained you in preparing and giving the injection.
- Ask your healthcare provider how often you will need to give the injection.
- Do not shake the pre-filled pen at any time.
- Do not remove the cap until you are ready to inject .
- Do not share the pre-filled pen with anyone.
- · Only use each pre-filled pen for one injection.

Prepare for the injection



GATHER THE SUPPLIES

- Prepare a clean, flat and well-lit surface such as a table.
- Remove the pen from the carton stored in your refrigerator.
- Items needed for injection:

Pen, Alcohol swab, Cotton ball or gauze, Adhesive bandage, Sharps disposal container

* Items not included in the carton

Give the injection



Medicine

Before use

Yuflyma

Plunger

After use

Pre-Filled Pen



- · Check the correct medicine and dosage.
- Check the expiration date on the label of the pen.
- colorless to pale brown and free of particles.
- * Do not use the pen if the liquid is discolored, cloudy, or contains particles.
- * It is normal for air bubbles to be seen in the liquid

WAIT 15 TO 30 MINUTES

- Leave the pen at room temperature (20°C to 25°C) for 15 to 30 minutes to let it warm up.
- * Do not warm the pen using heat sources.

INSPECT THE PEN AND THE MEDICINE · Look through the window and make sure the liquid is clear, * Do not use the pen if it has been cracked or is visibly damaged.

GIVE THE INJECTION • Press the pen firmly against the skin.

PLACE THE PEN

0 degree angle.

• When the injection starts you will hear the 1st loud and the blue plunger rod will begin to fill

the pen so that you can see the window.


thout pinching or stretching the skin,

e the pen over the injection site at

- Keep holing the pen and listen for the 2nd loud "click".
- · After the 2nd loud "click," continue to hold the pen and count slowly to 5.
 - * Do not change the position of the pen after the injection has started.



CHOOSE AN INJECTION SITE

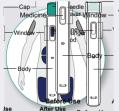

■ Caregiver ONLY

■ Self-injection

Select an injection site.

Possible injection sites:


- ✓ The front of the thighs
- ✓ Abdomen, except for the 5cm around the belly button(navel)
- ✓ Outer area of the upper arms, ONLY if given by a caregiver
- * Do not inject into skin that is within 5cm of your belly button (navel) or is tender, damaged, bruised, or scarred
- · Rotate the injection site each time you give an injection.
- Each new injection site should be at least 3cm away from the previous site.



REMOVE THE PEN

- Make sure that the blue plunger rod with the grey top has filled the window completely.
- Remove the pen from your skin.
- Gently press a cotton ball or gauze over the injection site.
- Do not reuse the pen.
- * Do not rub the injection site
- * If the window has not turned completely blue, this means you have not received a full dose. Call your healthcare provider immediately.

After the injection

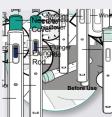

WASH YOUR HANDS AND CLEAN THE INJECTION SITE

- Wash your hands with soap and water and dry them fully.
- · Clean the injection site with an alcohol swab.
- Let the skin dry before injecting.
- * Do not blow or touch the injection site again before giving the injection.

DISPOSE OF THE PEN

- Throw away the used pen in a special sharps disposal container immediately after use.
- Alcohol pads and packaging may be put in your household waste.
- *Do not recap the pen.
- * Keep the pen and sharps disposal container away from children.

How do I inject with the Pre-filled Syringe?


Yuflyma is offered in an easy-to-use syringe. The following steps will guide you to a safe and easy administration of Yuflyma®.

Important information

- Use the pre-filled syringe only if your healthcare provider has trained you in preparing and giving the injection.
- Ask your healthcare provider how often you will need to give the injection.
- Do not shake the syringe at any time.
- Do not remove the cap until you are ready to inject .
- Do not share the syringe with anyone.
- · Only use each syringe for one injection.

Give the injection

REMOVE THE CAP

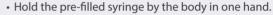

- · Pull the cap straight off.
- * Do not remove the cap until you are ready to inject.
- * Do not touch the needle.
- * Dispose of the cap immediately into the sharps disposal container.

Yuflyma

adalimumab

Pre-Filled Syringe

with Needle Guard

Finger

Flange

Window Medicine

After use

Before use

Prepare for the injection

GATHER THE SUPPLIES

- Prepare a clean, flat and well-lit surface such as a table.
- Remove the syringe from the carton stored in your refrigerator.
- Items needed for injection:

Pen, Alcohol swab, Cotton ball or gauze, Adhesive bandage,* Sharps disposal container*

* Items not included in the carton



INSPECT THE PEN AND THE MEDICINE

- Check the correct medicine and dosage.
- · Check the expiration date on the label of the syringe.
- Look through the window and make sure the liquid is clear, colorless to pale brown and free of particles.
- * Do not use the pen if it has been cracked or is visibly damaged.
- * Do not use the pen if the liquid is discolored, cloudy, or contains particles.
- * It is normal for air bubbles to be seen in the liquid

Step 2

Step 1



INSERT THE SYRINGE

- · Gently pinch a fold of skin at the injection site.
- While holding the syringe by its body, insert the needle completely into the fold of the skin at a 45-degree angle with a quick and "dart-like" motion.




WAIT 15 TO 30 MINUTES

- Leave the syringe at room temperature (20°C to 25°C) for 15 to 30 minutes to let it warm up.
- * Do not warm the pen using heat sources.

GIVE THE INJECTION

- After the needle is inserted, let go of the pinched skin.
- Push the plunger down slowly all the way down until the syringe is empty.
- * **Do not** change the position of the pen after the injection has started.

CHOOSE AN INJECTION SITE

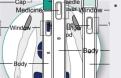

■ Caregiver ONLY


■ Self-injection


Select an injection site.

Possible injection sites:

- ✓ The front of the thighs
- ✓ Abdomen, except for the 5cm around the belly button(navel)
- ✓ Outer area of the upper arms, ONLY if given by a caregiver
- * Do not inject into skin that is within 5cm of your belly button (navel) or is tender, damaged, bruised, or scarred
- · Rotate the injection site each time you give an injection.
- Each new injection site should be at least 3cm away from the previous site.

REMOVE THE SYRINGE

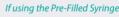


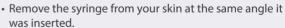

If using the Pre-Filled Syringe with Needle Guard:

- Slowly lift your thumb from the plunger until the needle is completely covered by the automatic needle guard.
- Gently press a cotton ball or gauze over the injection site and apply an adhesive bandage if necessary.
- * **Do not** rub the injection site.

After the injection

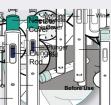

WASH YOUR HANDS AND CLEAN THE INJECTION SITE

- Wash your hands with soap and water and dry them fully.
- · Clean the injection site with an alcohol swab.
- Let the skin dry before injecting.
- * Do not blow or touch the injection site again before giving the injection.

DISPOSE OF THE SYRINGE

- Throw away the used syringe in a special sharps disposal container immediately after use.
- · Alcohol pads and packaging may be put in your household waste.
- * Do not recap the syringe.
- * Keep the syringe and sharps disposal container away from children.